

Phonak RemoteControl

User Guide



A Sonova brand


PHONAK
life is on

Welcome

Congratulations on choosing the Phonak RemoteControl.

Please read the user guide carefully to make sure that you understand and get the best out of your remote control.

For more information about features and benefits, simply contact your hearing care professional.

 Compatibility information:
Check your hearing aid compatibility with your hearing care professional.

Phonak – life is on
www.phonak.com




CE mark applied: 2019

Content

| | | | |
|---|----|--|----|
| 1. Description | 6 | 6. Compliance information | 16 |
| 2. Using Phonak RemoteControl | 7 | 7. Information and description of symbols | 22 |
| 2.1 Activating the Phonak RemoteControl | 7 | 8. Important safety information | 26 |
| 2.2 Switching On/Off | 7 | 8.1 Hazard warnings | 26 |
| 2.3 Understanding the indicator light | 7 | 8.2 Information on product safety | 29 |
| 2.4 Pairing the Phonak RemoteControl with your hearing aid | 8 | 8.3 Other important information | 31 |
| 2.5 Changing hearing aid volume | 9 | 9. Service and warranty | 32 |
| 2.6 Changing hearing aid program | 10 | 9.1 Local warranty | 32 |
| 3. Inserting a new battery | 11 | 9.2 International warranty | 32 |
| 3.1. Open the battery compartment | 11 | 9.3 Warranty limitation | 33 |
| 3.2. Remove the battery from the battery compartment | 11 | | |
| 3.3. Insert the new battery | 11 | | |
| 3.4. Close the battery compartment | 12 | | |
| 4. Reset the Phonak RemoteControl | 13 | | |
| 5. Troubleshooting | 14 | | |

1. Description


The intended user of the Phonak RemoteControl is to enable changes of your hearing aid programs and volume.

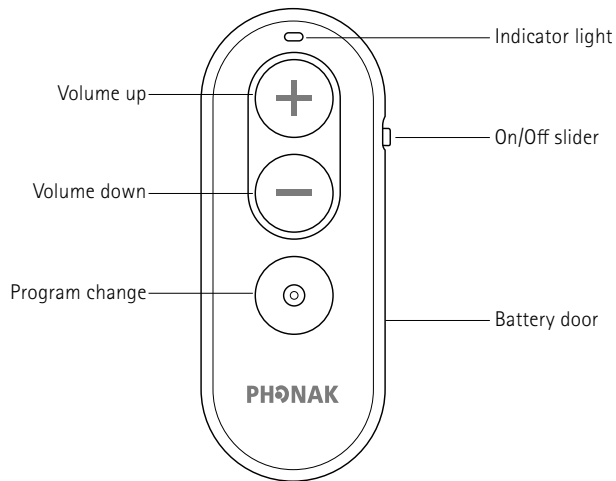

① The Phonak RemoteControl must be paired with your hearing aids before you can use it (chapter 2.4).


2. Using Phonak RemoteControl

2.1 Activating the Phonak RemoteControl
Pull out the battery protective tab to activate your Phonak RemoteControl.

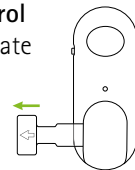





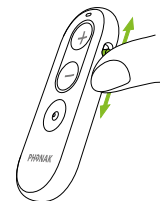



2.2 Switching On/Off

Use the slider to switch on or off.



2.3 Understanding the indicator light

| | | |
|---------------------------|---|-----------------------------------|
| Solid green for 3 seconds |  | Switched on, ready to use |
| Blinking blue |  | Pairing mode |
| Solid green for 5 seconds |  | Pairing successful |
| Solid red for 5 seconds |  | Pairing failed |
| Solid green for 1 second |  | Confirmation of push button press |
| Solid orange |  | Battery warning |
| Solid red |  | General device error |
| Solid green for 3 seconds |  | Reset successful |

2.4 Pairing the Phonak RemoteControl with your hearing aids


To use your Phonak RemoteControl with your hearing aids, you must first pair the hearing aids.

① It is only necessary to perform the pairing procedure once with your hearing aids. After the initial pairing, your Phonak RemoteControl will connect automatically with your hearing aids.


① The Phonak RemoteControl will be automatically in the pairing mode if switched on for the first time. This is indicated by the blinking blue light indicator on your Phonak RemoteControl.

Set your hearing aids into the pairing mode by switching them on. You now have 3 minutes to pair your Phonak RemoteControl with your hearing aids.

1. Press and hold the volume up "+" button
2. Switch on your Phonak RemoteControl whilst still pressing the volume up "+" button
3. Release the volume up "+" button when the light indicator starts to blink blue
4. The Phonak RemoteControl and your hearing aids will now pair automatically
5. After a successful pairing the light indicator is solid green for 5 seconds and you may hear a notification beep in your hearing aid
6. You can now use your Phonak RemoteControl

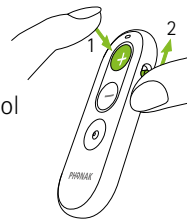


2.5 Changing hearing aid volume

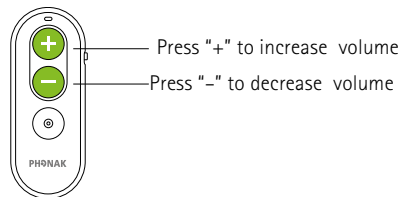

2.6 Changing hearing aid program

Each press of the program button will switch your hearing aids to the next available program.

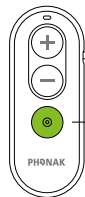

Press to change the hearing aid program


① Hearing aids may beep to confirm selection, depending on your hearing aid configuration.

① Changes are applied to both hearing aids at the same time.


① For information about your personal settings, ask your hearing care professional to print your Phonak RemoteControl instructions.

3. Inserting a new battery


3.1. Open the battery compartment

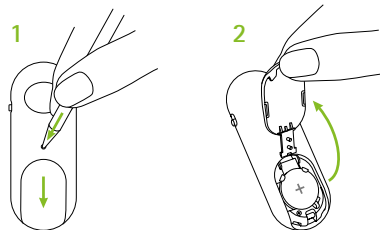

3.2. Remove the battery from the battery compartment

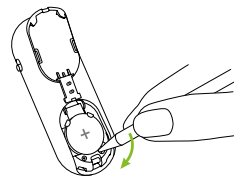


3.3. Insert the new battery

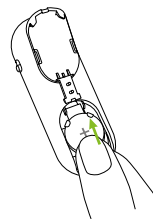

① The Phonak RemoteControl requires a lithium button cell (CR2032)

3.4. Close the battery compartment

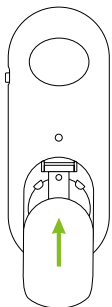

4. Reset the Phonak RemoteControl

① A reset of your Phonak RemoteControl will remove the pairing to your hearing aids.

Press the volume up "+" button and the program change button for more than 10 seconds.

5. Troubleshooting

Problem

No indicator light when a button is pressed

Hearing aids do not respond to the RemoteControl commands

Solid red light indicator

Causes

Remote control is switched off

Battery is empty

Hearing aids are switched off

Hearing aids are not paired with your RemoteControl

Remote control is switched off

Battery is empty

General device error

What to do

Switch on your remote control (chapter 2.2)

Change battery (chapter 3)

Switch on your hearing aids (see hearing aid user guide)

Pair the RemoteControl (chapter 2.4)

Switch on your remote control (chapter 2.2)

Change battery (chapter 3)

Reset your RemoteControl (chapter 4) and pair your hearing aids (chapter 2.4)

① If the problem persists, contact your hearing care professional for assistance.

6. Compliance information

Declaration of conformity

Hereby Phonak AG declares that this Phonak product is in compliance with the essential requirements of the Medical Devices Directive 93/42/EEC as well as the Radio Equipment Directive 2014/53/EU. The full text of the declaration of conformity can be obtained from the manufacturer or the local Phonak representative whose address can be taken from the list on <http://www.phonak.com/us/en/certificates> (worldwide locations).

Australia/New Zealand:



R-NZ

Indicates a device's compliance with applicable Radio Spectrum Management's (RSM) and Australian Communications and Media Authority (ACMA) regulatory arrangements for the legal sale in New Zealand and Australia. The compliance label R-NZ is for radio products supplied in the New Zealand market under conformity level A1.

The following model is certified under:

Phonak RemoteControl

USA FCC ID: KWC-RC1

Canada IC: 2262A-RC1

Notice 1:

This device complies with Part 15 of the FCC Rules and with RSS-210 of Industry Canada. Operation is subject to the following two conditions:

- 1) this device may not cause harmful interference, and
- 2) this device must accept any interference received, including interference that may cause undesired operation.

Notice 2:

Changes or modifications made to this RemoteControl not expressly approved by Phonak may void the FCC authorization to operate this device.

Notice 3:

This device has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules and ICES-003 of Industry Canada. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This device generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.


However, there is no guarantee that interference will not occur in a particular installation. If this device does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the device and receiver.
- Connect the device into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Notice 4:

Japanese Radio Law and Japanese Telecommunications Business Law Compliance

This device is granted pursuant to the Japanese Radio Law (電波法) and the Japanese Telecommunications Business Law (電気通信事業法). This device should not be modified (otherwise the granted designation number will become invalid).

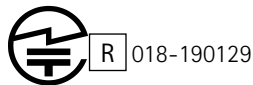

Radio information of your wireless device

| | |
|---------------------|--------------------|
| Antenna type | Monopole antenna |
| Operation frequency | 2.4 GHz – 2.48 GHz |
| Modulation | GFSK |
| Radiated power | < 2.8mW |
| Bluetooth® | |
| Range | ~1m |
| Bluetooth | 4.2 |
| Profiles supported | BLE (GATT) |

Compliance with emission and immunity standards

| | |
|--------------------|------------------------|
| Emission standards | EN 60601-1-2:2015 |
| | IEC 60601-1-2:2014 |
| | EN 55011:2009+A1 |
| | CISPR11:2009/AMD1:2010 |
| | CISPR22:1997 |
| | CISPR32:2012 |
| | ISO 7637-2:2011 |
| | CISPR25:2016 |
| | EN 55025:2017 |

| | |
|--------------------|-----------------------------|
| Immunity standards | EN 60601-1-2:2015 |
| | IEC 60601-1-2:2014 |
| | EN 61000-4-2:2009 |
| | IEC 61000-4-2:2008 |
| | EN 61000-4-3:2006+A1+A2 |
| | IEC 61000-4-3:2006+A1+A2 |
| | EN 61000-4-4:2012 |
| | IEC 61000-4-4:2012 |
| | EN 61000-4-5:2014 |
| | IEC 61000-4-5:2014 |
| | EN 61000-4-6:2014 |
| | IEC 61000-4-6:2013 |
| | EN 61000-4-8:2010 |
| | IEC 61000-4-8:2009 |
| | EN 61000-4-11:2004 |
| | IEC 61000-4-11:2004 |
| | IEC 60601-1 (§ 4.10.2):2005 |
| | ISO 7637-2:2011 |

7. Information and description of symbols



With the CE symbol, Sonova AG confirms that this product – including accessories – meets the requirements of the Medical Devices Directive 93/42/EEC as well as the Radio Equipment Directive 2014/53/EU. The numbers after the CE symbol correspond to the code of certified institutions that were consulted under the above-mentioned directives.



This symbol indicates that the products described in these user instructions adhere to the requirements for an applied part of Type B of EN 60601-1.



Indicates the medical device manufacturer, as defined in EU Directive 93/42/EEC.

Operating conditions

The product is designed such that it functions without problems or restrictions if used as intended, unless otherwise noted in these user guides. The remote control is classified IP52 and designed to be used in all daily life situations.



This symbol indicates that it is important for the user to read and take into account the relevant information in these user guides.



This symbol indicates that it is important for the user to pay attention to the relevant warning notices in these user guides.



Important information for handling and product safety.



This symbol confers that the electromagnetic interference from the device is under limits approved by the US Federal Communications Commission.



The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by Sonova AG is under license. Other trademarks and trade names are those of their respective owners.



Japanese mark for certified radio equipment.



Indicates the manufacturer's serial number so that a specific medical device can be identified.



Indicates the manufacturer's catalogue number so that the medical device can be identified.



This symbol indicates that it is important for the user to read and take into account the relevant information in this user guide.



Temperature during transportation and storage: -20° to $+60^{\circ}$ Celsius (-4° to $+140^{\circ}$ Fahrenheit).



Humidity during storage: 0% to 70%, if not in use.



Atmospheric pressure during transportation and storage: 200 hPa to 1500 hPa



During transportation keep dry.








The symbol with the crossed-out garbage bin is to make you aware that this device may not be thrown away as normal household waste. Please dispose of old or unused device, at waste disposal sites intended for electronic waste, or give your device to your hearing care professional for disposal. Proper disposal protects the environment and health.


8. Important safety information


Please read the information on the following pages before using your device


8.1 Hazard warnings





-  Your device operate between 2.4 GHz–2.48 GHz frequency range. When flying please check if flight operator requires devices to be switched off.
-  The intended use of the remote control is to enable changes of the hearing aid programs and volume.
-  Keep this device out of reach of children and mentally challenged people or pets.
-  This remote control uses low-power, digitally coded transmission to communicate to your hearing system. Although unlikely, interference with medical devices (i.e. pacemakers, defibrillators, etc.) is possible. Therefore users must not keep this remote control in or near a breast pocket.

 Changes or modifications to the device that were not explicitly approved by Sonova are not permitted. Such changes may damage your device.





 If problems occur which cannot be resolved by following the remedy guidelines in the troubleshooting section of this user guide, consult your hearing care professional.

 Batteries are toxic if they are swallowed! Keep out of the reach of children, individuals with cognitive impairment, and pets. If batteries are swallowed, consult your physician immediately!

 This device is not for children below 36 months. The usage of this device by children and individuals with cognitive impairment should be supervised at all times to ensure their safety. It is a small device and contains small parts. Do not leave children and individuals with cognitive impairment unsupervised with this device. If swallowed, consult a physician or hospital immediately as the device or its parts can cause choking!

-  Remove the battery if you are not using the device for a long period of time.
-  Dispose of electrical components in accordance with your local regulations by Sonova AG.
-  External devices may only be connected if they have been tested in accordance with corresponding IECXXXXX standards.
-  Do not use the device in explosive areas (flammable anesthetics, mines or industrial areas with a danger of explosions), oxygen rich environments or where electronic equipment is prohibited.

8.2 Information on product safety

-  Protect the device from excessive moisture (bathing, swimming), heat (radiator, car dashboard) or direct skin contact when sweating (workout, fitness, sport).
-  Do not drop the device. Dropping onto a hard surface can damage your device.
-  Special medical or dental examination including radiation described below, may adversely affect the correct functioning of your device. Remove and keep it outside the examination room/area before undergoing:
 - Medical or dental examination with X-ray (also CT scan).
 - Medical examinations with MRI/NMRI scans, generating magnetic fields.
-  Protect the device sockets and battery slot from dirt and debris.

- ① Clean the device using a damp cloth. Never use household cleaning products (washing powder, soap, etc.) or alcohol to clean the device. Never use a microwave or other heating devices to dry the device.
- ① The digitally-coded, transmission technology used in the device is highly reliable and experiences virtually no interference from other devices. It should be noted, however, that when operating the hearing system near computer equipment, larger electronic installations or other strong electromagnetic fields, it may be necessary to be at least 60 cm (24") away from the interfering device to ensure proper operation.
- ① Keep the device at least 10 cm away from any kind of magnets.
- ① When the device is not in use, turn it off and store it safely.
- ① Do not use your device in areas where electronic equipment is prohibited.

8.3 Other important information

- ① High-powered electronic equipment, larger electronic installations and metallic structures may impair and significantly reduce the operating range.
- ① If the hearing aids do not respond to the device because of an unusual field disturbance, move away from the disturbing field.

9. Service and warranty

9.1 Local warranty

Please ask the hearing care professional, where you purchased your devices, about the terms of the local warranty.

9.2 International warranty

Sonova offers a one year limited international warranty, valid as of the date of purchase. This limited warranty covers manufacturing and material defects. The warranty is valid only if proof of purchase is shown.

The international warranty does not affect any legal rights that you might have under the local warranty or the applicable national legislation governing sale of consumer goods.

9.3 Warranty limitation

This warranty does not cover damage from improper handling or care, exposure to chemicals or undue stress. Damage caused by third parties or non-authorized service centers renders the warranty null and void. This warranty does not include any services performed by a hearing care professional in their office.

Serial number: _____

Date of purchase: _____

Authorized hearing care professional
(stamp/signature):

Your hearing care professional:



Manufacturer:

Sonova AG
Laubisrütistrasse 28
CH-8712 Stäfa
Switzerland
www.phonak.com



029-0782-02/V1.00/2019-08/na © 2019 Sonova AG All rights reserved

sonova
HEAR THE WORLD

